


GEISTER

GEISTER®

The better way to operate™

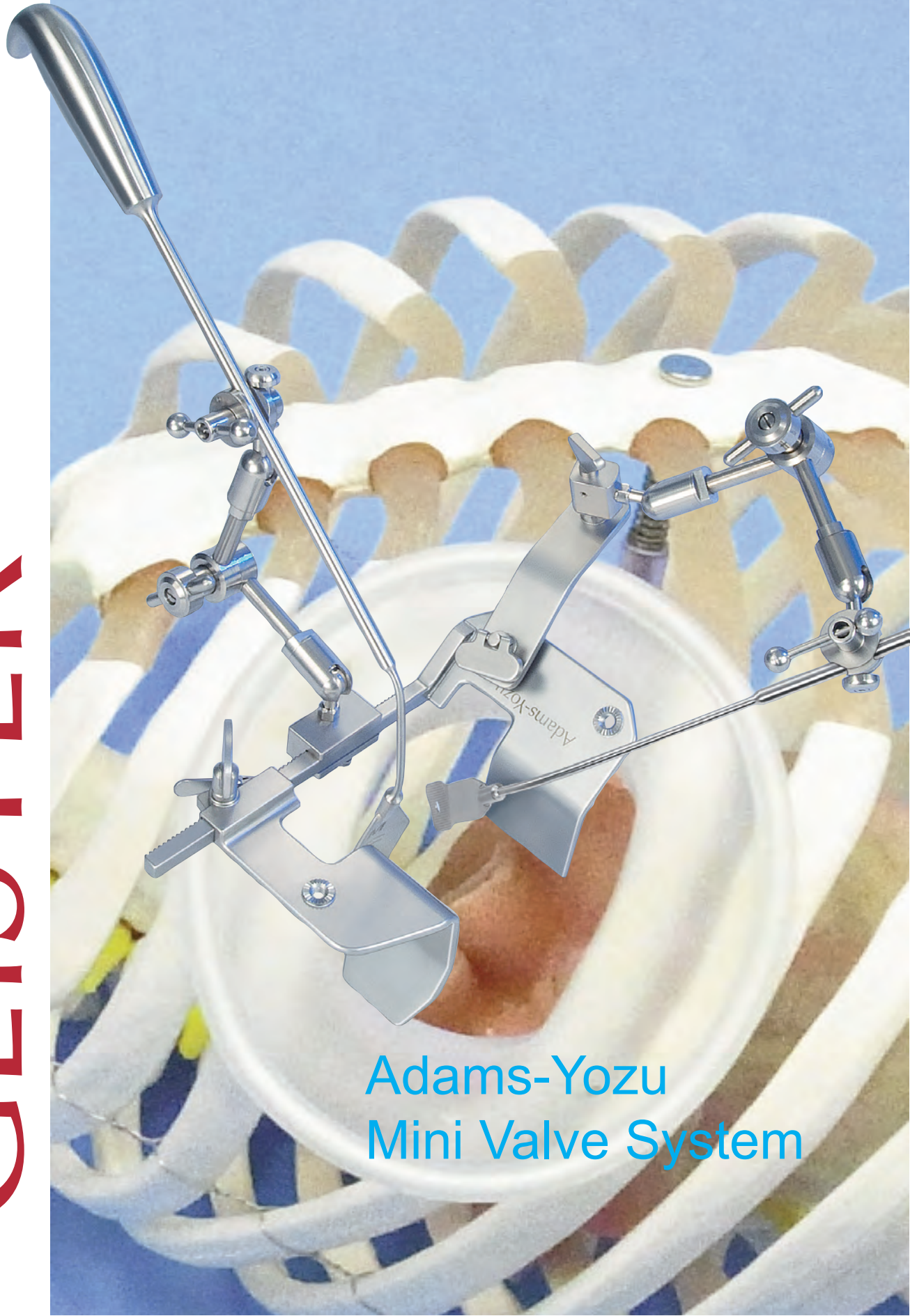

Adams-Yozu
Mini Valve System

This retractor system facilitates both direct-vision and endoscopic-assisted mini-mitral valve surgery and helps the surgeon to accomplish a favorable surgical outcome even in the complex pathology of mitral valve. Special left atrial blades and an additional blade for the posterior wall of the left atrium can be attached to the Adams-Yozu retractor. The blades can be flexibly moved and then fixed in any favorable position to achieve optimal exposure of the mitral valve. 5 blade sizes are available depending on the left atrial size and target lesion. By using the smallest size, even papillary muscles can be exposed easily and clearly. The angled chest spreader creates wider working space around the ascending aorta with opening. This helps surgeon's manipulation of cross-clamping and the aortic root.



29-0930




29-0934



29-0933



29-0935

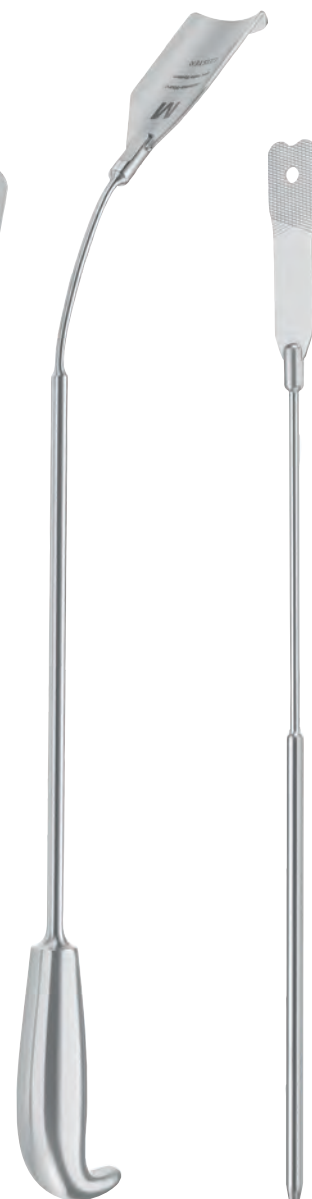

29-0937.35

29-0938

Order information:

29-0930 Adams-Yozu Minithoracotomy Retractor

Accessories:

- 29-0933 Adams-Yozu Mini Arm with M6 screw connector
- 29-0934 Adams-Yozu Mini Arm with quick connector for spreading bar
- 29-0935 Attachment and Extension blade
- 29-0937.23 Depressor blade 23mm „XS“, 70mm deep
- 29-0937.28 Depressor blade 28mm „S“, 70mm deep
- 29-0937.35 Depressor blade 35mm „M“, 70mm deep
- 29-0937.45 Depressor blade 45mm „L“, 70mm deep
- 29-0937.45L Depressor blade 45mm „XL“, 90mm extra deep
- 29-0938 Malleable Spatula



The better way to operate™

Geister Medizintechnik GmbH
 Foehrenstraße 2
 D-78532 Tuttlingen/Germany
 Phone +49-(0)74 61-96 62 40
 Telefax +49-(0)74 61-96 62 422
 E-mail: info@geister.com
 Internet: www.geister.com

Distributor